


# TECHNICAL DATA

DESCRIPTION	351 SHORT FRAME	352 MEDIUM FRAME	353 LONG FRAME
Seat width	300 - 360	360 - 580	360 - 580
Width when ready for use	Seat width + 200 (with side guard CODE 100 SB + 180)		
Folded width	280		
Seat depth	330 - 380	350 - 460	480 - 530
Seat height front	370 - 500	370 - 560	370 - 560
Seat height back	400 - 470	390 - 520	390 - 520
Overall height	830	920	920
Back height	340 - 400	340 - 500	340 - 500
Armrest height (CODE 81)	165 - 240	200 - 290	200 - 290
Length with footrests	880	1040	1120
User weight	75	135	160
Min. chair weight	12	13	14
Permitted total weight	87	148	174
Transport weight	8	8.5	9
Front wheel (22")	560	560	560
Front wheel (24")		610	610

Dimensions in mm, weight in kg, wheels in inches, speed in km/h, size tolerance  $\pm 10$  mm, data subject to changes in design.

## FRAME COLOURS

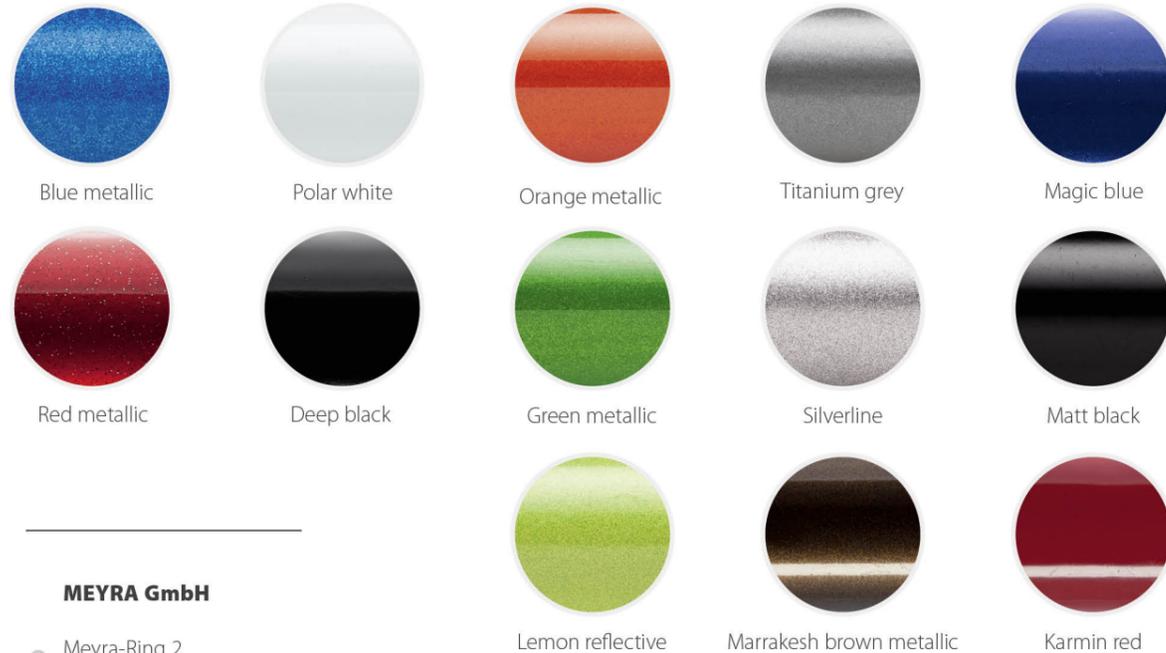

### MEYRA GmbH

Meyra-Ring 2  
D-32689 Kalletal-Kalldorf

Phone +49 5733 922 0  
Fax +49 5733 922 9311

info@meyra.de

www.meyra.de

by www.bwc.gmbh



MEYRA®

**AVANTI**

OPTIMISED TO THE LAST DETAIL

ACTIVE

# » DISCOVER THE NEW POSSIBILITIES «

ALSO NEW

- 13 frame colours at no extra cost for customized colour design, Vario block and castor fork available in different frame colour on request
- Back strap height available in 460, 480 and 500 mm
- Seat belt can be tightened
- 50 mm wheel cover
- Aluminium double cross brace for seat widths between 360 - 460 mm, recommended for use with auxiliary drives or under high stress
- Swing-away or mountable anti-tip wheels, can be used as a pair or individually



Seat width in 420 mm and 440 mm



Also available with high-adjustable side guard, can be operated with one hand



Telescopic and angle-adjustable back height from 350 to 410 mm as well as 400 to 460 mm for ideally adapting to changing symptoms



Compact (78° anterior frame angle) or standard (70° anterior frame angle) upper part of leg rest, removable and swing-away



Lower leg length can be adjusted in steps or progressively to meet the user's requirements, a shorter lower leg length of 290 - 390 mm is also available



Up to 12° seat inclination for precise adjustment options



Active wheel with radial spokes, Schwalbe ONE tyres, hub, rim, hand rims in blue, red or yellow for more individualised options



Schwalbe ONE tyres, no profile, puncture-resistant and lightweight



Special brake with extremely low-resistance activation



Curve GEKKO, black anodised aluminium with replaceable rubber lip, ergonomic grip profile, ideal for tetraplegic users